



Quantitative measurement and visualization of biofilm O₂ consumption rates inmembrane filtration systems

Vrouwenvelder, J. S.; Prest, E. I.E.C.; Kuhl, M.; Van Loosdrecht, M. C.M.; Staal, M.

Published in:
Procedia Engineering

DOI:
[10.1016/j.proeng.2012.08.370](https://doi.org/10.1016/j.proeng.2012.08.370)

Publication date:
2012

Document version
Publisher's PDF, also known as Version of record

Document license:
[CC BY-NC-ND](#)

Citation for published version (APA):
Vrouwenvelder, J. S., Prest, E. I. E. C., Kuhl, M., Van Loosdrecht, M. C. M., & Staal, M. (2012). Quantitative measurement and visualization of biofilm O₂ consumption rates inmembrane filtration systems. *Procedia Engineering*, 44, 233-234. <https://doi.org/10.1016/j.proeng.2012.08.370>

Euromembrane Conference 2012
[OB20]
**Quantitative measurement and visualization of biofilm O₂ consumption rates
inmembrane filtration systems**

 J.S. Vrouwenvelder^{*1,2}, E.I.E.C. Prest^{1,3}, M. Kuhl^{4,5}, M.C.M. van Loosdrecht¹, M. Staal^{2,4}
¹*Delft University of Technology, The Netherlands*, ²*KAUST, Saudi Arabia*, ³*Wetsus, The Netherlands*, ⁴*University of Copenhagen, Denmark*, ⁵*University of Technology Sydney, Australia*

There is a strong need for techniques enabling direct assessment of biological activity of biofouling in membrane filtration systems. Here we present a new quantitative and non-destructive method for mapping O₂ dynamics in biofilms during biofouling studies in membrane fouling simulators (MFS). Transparent planar O₂ optodes in combination with a luminescence lifetime imaging system were used to map the two-dimensional distribution of O₂ concentrations and consumption rates inside the MFS. The O₂ distribution was indicative for biofilm development. Biofilm activity was characterized by imaging of O₂ consumption rates, where low and high activity areas could be clearly distinguished. The spatial development of O₂ consumption rates, flow channels and stagnant areas could be determined (Figures 1 and 2). This can be used for studies on concentration polarization, i.e. salt accumulation at the membrane surface resulting in increased salt passage and reduced water flux. The new optode-based O₂ imaging technique applied to MFS allows non-destructive and spatially resolved quantitative biological activity measurements for on-site biofouling diagnosis and laboratory studies. The following set of complementary tools is now available to study development and control of biofouling in membrane systems: (i) MFS, (ii) sensitive pressure drop measurement, (iii) magnetic resonance imaging, (iv) numerical modelling, and (v) biological activity measurement based on O₂ imaging methodology.

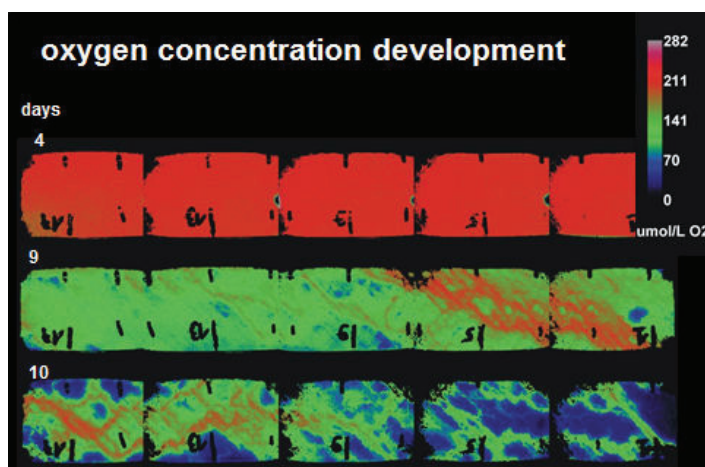


Figure 1: Spatio-temporal distribution of O₂ concentration (μmol O₂ / L) imaged over the monitor length at the optode surface, i.e. the base of the fouling layer. The figure shows the development of low O₂ concentration regions and flow channelling over time. The flow direction is from left to right.

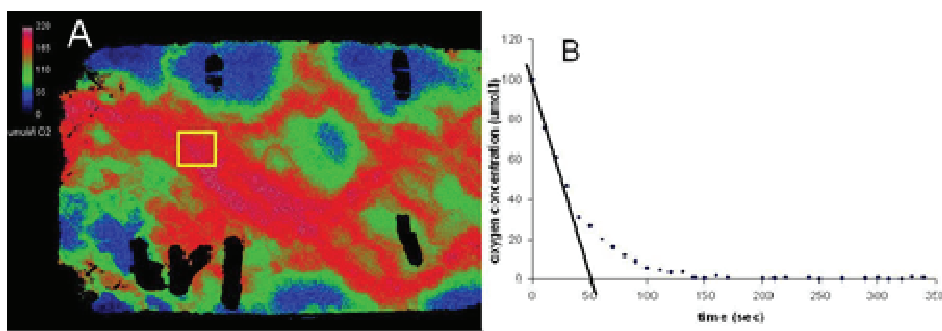


Figure 2. Procedure to determine O_2 consumption rates in the MFS. Oxygen images were taken every 10 s after the flow was stopped. An area of interest (AOI) was selected on the initial oxygen image (A) and the O_2 depletion over time was determined for the AOI (B). The O_2 consumption rate for particular AOIs was calculated from the initial slope of the O_2 depletion curve.

Keywords: biofouling, non-destructive biofouling diagnosis, biological activity measurement, concentration polarization